

DTIC FILE COPY

4

AD

AD-E402 089

Contractor Report ARAED-CR-90011

AD-A226 186

**DEVELOPMENT OF AQUEOUS AND NONAQUEOUS  
TITRIMETRIC METHODS FOR LIQUID PROPELLANTS**

J.P. Griff  
G. Doyle  
Geo-Centers Inc.  
762 Route 15 South  
Lake Hopatcong, NJ 07849

DTIC  
ELECTE  
SEP 06 1990  
S D 5 D

E. Turngren  
W.O. Seals  
Project Engineers  
ARDEC

August 1990



US ARMY  
ARMAMENT MUNITIONS  
& CHEMICAL COMMAND  
ARMAMENT RDE CENTER

**U.S. ARMY ARMAMENT RESEARCH, DEVELOPMENT AND  
ENGINEERING CENTER**

Armament Engineering Directorate  
Picatinny Arsenal, New Jersey

Approved for public release; distribution is unlimited.

90 09 05 062

The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other documentation.

The citation in this report of the names of commercial firms or commercially available products or services does not constitute official endorsement by or approval of the U.S. Government.

Destroy this report when no longer needed by any method that will prevent disclosure of contents or reconstruction of the document. Do not return to the originator.

UNCLASSIFIED  
SECURITY CLASSIFICATION OF THIS PAGE

## REPORT DOCUMENTATION PAGE

1a. REPORT SECURITY CLASSIFICATION <b>UNCLASSIFIED</b>		1b. RESTRICTIVE MARKINGS													
2a. SECURITY CLASSIFICATION AUTHORITY		3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for public release; distribution is unlimited													
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE															
4. PERFORMING ORGANIZATION REPORT NUMBER		5. MONITORING ORGANIZATION REPORT NUMBER Contractor Report ARAED-CR-90011													
6a. NAME OF PERFORMING ORGANIZATION Geo-Centers	6b. OFFICE SYMBOL	7a. NAME OF MONITORING ORGANIZATION ARDEC, AED													
6c. ADDRESS (CITY, STATE, AND ZIP CODE) 762 Route 15 South Lake Hopatcong, NJ		7b. ADDRESS (CITY, STATE, AND ZIP CODE) Energetics and Warheads Division (SMCAR-AEE) Picatinny Arsenal, NJ 070806-5000													
8a. NAME OF FUNDING/SPONSORING ORGANIZATION ARDEC, IMD STINFO Br	8b. OFFICE SYMBOL SMCAR-IMI-I	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER DAAA 21-88-D-0021													
8c. ADDRESS (CITY, STATE, AND ZIP CODE) Picatinny Arsenal, New Jersey 070806-5000		10. SOURCE OF FUNDING NUMBERS <table border="1"><tr><td>PROGRAM ELEMENT NO.</td><td>PROJECT NO.</td><td>TASK NO.</td><td>WORK UNIT ACCESSION NO.</td></tr><tr><td></td><td></td><td></td><td></td></tr></table>		PROGRAM ELEMENT NO.	PROJECT NO.	TASK NO.	WORK UNIT ACCESSION NO.								
PROGRAM ELEMENT NO.	PROJECT NO.	TASK NO.	WORK UNIT ACCESSION NO.												
11. TITLE (INCLUDE SECURITY CLASSIFICATION) Development of Aqueous and Nonaqueous Titrimetric Methods for Liquid Propellants															
12. PERSONAL AUTHOR(S) J.P. Griff and G. Doyle, Geo-Centers Inc., E. Turngren and W.O. Seals, ARDEC Project Engineers															
13a. TYPE OF REPORT Progress	13b. TIME COVERED FROM _____ TO _____	14. DATE OF REPORT (YEAR, MONTH, DAY) August, 1990	15. PAGE COUNT 34												
16. SUPPLEMENTARY NOTATION															
17. COSATI CODES <table border="1"><tr><th>FIELD</th><th>GROUP</th><th>SUB-GROUP</th></tr><tr><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td></tr></table>		FIELD	GROUP	SUB-GROUP										18. SUBJECT TERMS (CONTINUE ON REVERSE IF NECESSARY AND IDENTIFY BY BLOCK NUMBER) Aqueous    Liquid propellants    Nitric acid    Nonaqueous    Iron Hydroxylammonium nitrate    Surveillance    Triethanolammonium nitrate	
FIELD	GROUP	SUB-GROUP													
19. ABSTRACT (CONTINUE ON REVERSE IF NECESSARY AND IDENTIFY BY BLOCK NUMBER) <p>The U.S. Army has made a commitment to develop liquid gun propellants as a potential replacement of solid propellants used in various gun systems. In order to monitor the propellants during storage, a study was conducted to develop analytical capabilities for this purpose. The storability of liquid propellants is required to determine whether there is any instability, its cause, and if necessary, its prevention. Titrimetry was selected to monitor two of the species present which were most illustrative of propellant conditions and critical to performance. In the course of this study, the merits of aqueous and nonaqueous techniques were discovered. As a result of this study, a capability is available to monitor some of the vital components and contaminants in liquid propellants during storage.</p>															
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT <input type="checkbox"/> UNCLASSIFIED/UNLIMITED <input checked="" type="checkbox"/> SAME AS RPT. <input type="checkbox"/> DTIC USERS		21. ABSTRACT SECURITY CLASSIFICATION UNCLASSIFIED													
22a. NAME OF RESPONSIBLE INDIVIDUAL I. HAZNEDARI		22b. TELEPHONE (INCLUDE AREA CODE) (201) 724-2316	22c. OFFICE SYMBOL SMCAR-IMI-I												

# CONTENTS

	Page
Introduction	1
Experimental	2
Chemicals	2
Titrators	2
Results	2
Conclusion	8
References	19
Distribution List	21



Accession For	
NTIS CR441	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution /	
Availability Codes	
Dist	Avail and/or Special
A-1	

## TABLES

1 Determination of nitric acid spiked in propellant samples by nonaqueous titration	9
2 Determination of nitric acid in spiked propellant samples by aqueous titration	10
3 Effect of diluent blank in nitric acid study	11
4 Comparison of aqueous versus nonaqueous titration for HAN and TEAN	11
5 Comparison of aqueous and nonaqueous synthetic 1846 HACL/TEACL titrations	12

## FIGURES

	Page
1 Nitric acid titration % recovery versus sample weight	13
2 Nonaqueous titration derivative curves with optimized sample size	14
3 Nonaqueous titration derivative curve with concentrated sample	14
4 Effect of sample concentration in aqueous titration derivative curves	15
5 Aqueous titration derivative curves	16
6 Potentiometric titrations S-curve	16
7 Aqueous titration S-curve for nitric acid	17
8 Nonaqueous titrations S-curves for notric acid	17
9 Aqueous titration S-curve for HAN/TEAN	18
10 Nonaqueous titration S-curve HAN/TEAN	18

## INTRODUCTION

The U.S. Army has been engaged in the search for intensive munitions for guns in tanks and self-propelled howitzers. Liquid Propellants have been categorized as an intensive munition (ref 1) and are presently under consideration as a potential replacement of solid propellants. The logistical, ballistic, and cost savings feature propellant prompted an extensive program for their characterizations. The integrity of these propellants, after being stored for extended periods of time, is required to assure that the ballistic parameters are met. In order to establish these criteria, long-term storage studies of liquid propellants are required to determine their stability and storability.

The use of liquid propellants as a replacement for solid propellants in diverse gun applications necessitates extensive characterization of the propellant system. This requirement has resulted in the necessity to develop analytical techniques that will establish the capability to monitor the composition of liquid propellants during storage. Some of the important features which were considered for candidate techniques included reproducibility, reliability, simplicity, safety, and environmental impact.

The propellant systems currently under investigation are stoichiometric mixtures of hydroxylammonium nitrate (HAN) with an alkyl ammonium nitrate in water. The alkylammonium nitrate which is being considered at present is triethanolammonium nitrate (TEAN). A typical composition called LP 1846 has 61% HAN, 19% TEAN, and 20% water. The molar ratio of this mixture is 7:1 for HAN to TEAN which is stoichiometric for conversion to the desirable products of combustion of carbon dioxide, water, and nitrogen.

The knowledge of the aging effects of the basic liquid propellant constituents and stability when exposed to contaminants during long-term storage is essential. This information is critical since any significant deviation in the basic propellant formulation will result in a loss in ballistic performance and propellant lifetime as well as causing potential safety problems. Therefore, it is necessary to determine the effects of temperature and contaminants on the stability of liquid propellants and provide analytical techniques to monitor the gas and liquid phase compositions before and after exposure.

Analyses are required to both identify and quantify major propellants, contaminants, and degradation products in order to provide a basis for desired criteria. Low level concentrations of contaminants or impurities which have been detailed previously (ref 2 and 3) such as ammonium nitrate (AN), nitric acid, morpholines, hydrazines, nitrosoamines, nitrogen oxides, nitrogen, and trace amounts of transition metal ions can be expected to be found.

## EXPERIMENTAL

### Chemicals

LP 1846 propellant was prepared by Thiokol Chemical Corp., Elkton, MD and supplied by BRL, Aberdeen, MD.

The hydrochloride derivatives of HAN and TEAN were purchased from Aldrich Chemical Co., Milwaukee, WI.

The high purity solvents and reagents are commercial products. Ethanol (Pharmco, NJ), sodium hydroxide, methanol, and acetone (Fisher, NJ).

The tetrabutylammonium hydroxide is a commercially available titrant in methanol (Aldrich Chem Co.).


### Titration

The potentiometric titrations were carried out on a Metrohm potentiograph Model E-536 and a Model 655 Metrohm Dosimat.


## RESULTS

Several of the analyses have received considerable attention due to their complexity and necessity for accurate results. The determination of both free nitric acid and of hydroxylammonium nitrate (HAN) and triethanolammonium nitrate (TEAN) is an area that has provided some difficulty. In-depth investigations of the use of titrimetry as well as alternative techniques such as ion chromatography (IC) and supercritical fluid chromatography (SFC) have been conducted.

A titrimetric method developed by Dr. Kasler at the University of Maryland was reviewed as a potential technique for HAN/TEAN analysis (refs 4 and 5). Methyl isobutyl ketone was used in this method to convert HAN to a stronger acidic product, nitric acid. TEAN, a weak acid, is not affected by the ketone and is easily differentiated in the titration curve. The actual titration is based on the nitric acid liberated during the formation of an oxime, a very weak acid ( $K_b < 10^{-12}$ ) which is not titrated.

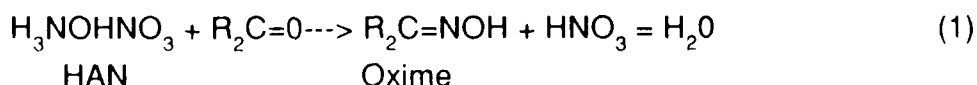

HAN based liquid propellants can undergo decomposition as a result of the presence of various contaminants such as nitric acid and transition metals with formation of numerous products some of which are shown in equation 2:

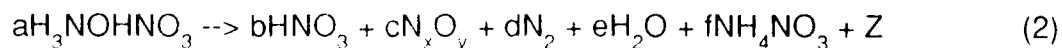

The actual mechanism for liquid propellant decomposition proposed in equation 2 will be established by gas and liquid phase analyses of experimentally aged samples. Nitric acid formation during decomposition poses a serious problem in the stability of the liquid propellant and must be accurately monitored. "Free" nitric acid can be introduced as a product of decomposition or from production as a residual contaminant. Therefore, this titrimetric method will result in erroneously high HAN values in samples with significant decomposition or with residual acid from production. As a result, a capability was required to differentiate free nitric acid separate from nitric acid produced from keytone/HAN reaction. This was accomplished by performing two titrations. The keytone was eliminated from the sample solvent permitting the HAN to remain unchanged. The free nitric acid from HAN to remain unchanged. The free nitric acid from HAN decomposition or propellant production could then be titrated. This is impossible since the HAN and TEAN which are very weak acids from one break in the titration curve after the stronger, free nitric acid break. Tetrabutylammonium hydroxide of 0.01N in methanol was used as the titrant with 50 ml of ethanol as the sample solvent. In order to provide the necessary accuracy for concentrations of  $\text{HNO}_3 < 0.5\%$ , an 0.6-gram sample was chosen; and for concentrations of  $\text{HNO}_3 > 0.5\%$ , an 0.3-gram sample was required (refs 2 and 3). Finally, the HAN/TEAN could be determined in a separate titration with addition of keytone to the sample solvent for HAN conversion.

A comparison of the aqueous versus nonaqueous free nitric acid titrations led to the selection of the nonaqueous method. The overall standard deviation of the aqueous was  $> +0.02$  and for the nonaqueous  $< +0.02$  (tables 1 and 2). The titration curves for the nonaqueous method were more sensitive. The choice of sample weight was also critical. In keeping with acceptable titration techniques as prescribed by numerous authors as well as safety and environmental directives, sample size was kept minimal. The effect of sample concentration is vividly expressed in figure 1 where percent recovery is dramatically reduced as sample concentration is increased. The calculations for figure 1 were made without subtraction of the solvent blank in order to show the effect of large samples. The actual recovery of nitric acid added (table 1) is a complete loss or zero recovery with sample sizes greater than one gram.

The choice of methods for this study was restricted by sample quantities available, safety regulations, directives to reduce waste and sample usage, and analytical methods had to be tailored to comply with these restrictions.



The original HAN/TEAN method by Kasler produced titration curves which produced breaks that were not sufficiently discernible for the accuracy required for this program. In this study, first derivative curves were used to provide more easily defined end points. Acetone was substituted for the higher keytones used by other investigators (refs 3 and 4) for conversion of HAN to nitric acid which led to smoother and better defined titration curves. Pre-addition of titrant was used to reduce analysis time and provide accurate data. The sample size (0.6 g) was within accepted levels for ionic strengths of less than 0.1 (eqs 5 through 7), but sufficient for the desired accuracy and precision. Several titrants were tested, and tetrabutyl ammonium hydroxide (TBAH) was found to be applicable to both nitric acid and HAN/TEAN analyses. The preferred titrant for free nitric acid from degradation and production was n-butylamine in methanol, but the use of both titrants was not pursued. It was expedient to use only one titrant since both gave comparable results. The titrant for the HAN/TEAN analytical method is 0.2N TBAH in methanol. The sample size is 0.6 g with pre-addition of 15.0 ml of titrant in 50 ml of ethanol/acetone (100:1) (refs 2 and 3).

Blank values were determined for all reagents. The effect of the reagent used as a diluent or solvent was in some cases quite significant. A tabulation of these data are shown in table 3. Calculation of the blank as nitric acid shows that a significant error would be introduced in the case of ethanol 2.

The aqueous system was compared to nonaqueous titration of two liquid propellant lots (table 4 ) and of synthetic lot using hydrochloride standards (table 5). Both LP-2 and LP-3 were specification grade mixtures containing  $60.8 \pm 0.5\%$  HAN. As a result, the data in table 4 indicate that the nonaqueous titration provides more realistic results. This is also confirmed with the data in table 5. In this case, pure hydrochloride standards were used to prepare solutions containing the same amount of cation (hydroxylamine or triethanolamine) which would be present in the nitrate solutions of LP 1846. Again, the nonaqueous determination resulted in closer agreement with actual concentration of HAN. The only interesting observation in the TEAN comparison is that the aqueous titration results in higher TEAN (TEAC1) concentrations. This may be due to unreacted HAN (HAC1) which would account for low HAN values. If the effect of molecular weight from equation 3 is considered, any unreacted HAN (MW=96) would have more than a two-fold effect on increasing the percent TEAN (MW=212) detected. Considering the role of molecular weight in the titrimetry calculations, the presence of a species of lower molecular weight would increase the calculated amount detected (eq 3). For example, a weight of 80 versus TEAN with a molecular weight of 212 would effectively increase the TEAN concentration detected due to its presence by more than twice the amount.

$$\%X = \frac{(100) [(mls \text{ titrant}) (N \text{ titrant}) (MW/1000)]}{\text{Sample weight, grams}} \quad (3)$$

In order to provide an analytical capability, an overall review and understanding of the basic principals of the methods employed is necessary. One of the approaches, titrimetry, which has been chosen for characterizing the major components and the nitric acid contaminant in liquid propellants is one of complex dimensions. Titrimetry in either aqueous or nonaqueous media is often neither simple nor well defined and must be used with an understanding of its limitations. For this particular application, nonaqueous potentiometric titrimetry was selected for reasons which will be enumerated.

The selection of nonaqueous over aqueous titration methods provides increased sensitivity of weak acidic species such as TEAN. That is, the amplitude of the derivative curve is greater and much sharper. By careful observance of titrant and solvent effects and the increased sensitivity of nonaqueous media, it was possible to qualitatively and semiquantitatively distinguish between weakly acidic species, TEAN and AN. This characteristic has been observed by others in ethanol and other nonaqueous media (ref 6).

To be consistent with general practices set forth by numerous investigators, a dilute system was considered as an optimum approach with ionic strengths less than 0.1 (ref 7). These criteria are fulfilled by both nonaqueous methods which were developed (eqs 5 through 7). Concentrated solutions create large changes in activity coefficients and buffering affecting the shape of the titration curves (refs 6 and 8) which was evidenced in recent studies of percent recovery of nitric acid with sample size (fig. 1). Since buffering capacity is the property of a solution to resist pH change, it is requisite to reduce buffering power through dilution especially in analyses which require monitoring traces or slight changes. Evidence of buffering as a result of concentrated samples is illustrated in the derivative curves of figures 2 through 5. A standard symmetrical nonaqueous titration is shown in figure 2; the effect of a larger sample in which the curve becomes very noisy and shallow in figure 3; and the effect of increasing sample size using the aqueous method whose curves are asymmetrical in figure 4. As sample size was increased, the curve became increasingly shallow and eventually indiscernible. Figure 5 is an aqueous titration with a very large sample (~30 grams) in which the curve is measurable but very noisy and asymmetrical. The agreement obtained between the aqueous and the nonaqueous methods may have been coincidental. It was a deviation from the other concentrated samples which were tried since a discernible break was observed. Nevertheless, the more reasonable approach is with small samples.

Using classical associations and definitions, the following relationship exists between pH and pK values:

$$\text{pH} = \text{pK}_a + \log (\text{A}_b/\text{A}_a) \text{ if } \text{A}_a = \text{A}_b, \text{pK}_a = \text{pH}_k \quad (4)$$

where  $\text{A}_a$  and  $\text{A}_b$  are activities ( $\sim$  concentrations).

Since  $\text{A}_a = \text{A}_b$  is approximately fulfilled by  $V_E/2$ , the point K on the curve can be easily determined resulting in the appropriate pH<sub>k</sub> form which the corresponding pK<sub>a</sub> can be

found (fig. 6). Although the pKa determined in this manner is not exact, it is sufficient for this comparison. The buffering capacity,  $b$ , is a differential quantity and can be determined by drawing a tangent to the titration curve (fig. 7). The buffering capacity can also be calculated as  $b = DV \text{ (equiv/L) } / Dph$ . The more symmetrical the curve, the more distant  $b$  is from the end point. Titrations were run using titration techniques currently being employed for HAN-based propellants (figs. 7 through 10). Calculations are shown below for several of the aqueous and nonaqueous titration S-curves of LP 1846-03-11 from figures 7 through 10:

Method	End point		$b$	pK	%HNO <sub>3</sub>	%HAN	Start
	ml	pH					pH
Aqueous	0.50	2.08	0.28	1.92	0.03		1.80
Nonaqueous	1.50	2.17	0.019	1.79	0.07		1.54
Aqueous	15.53	5.04	0.24	1.90		60.29	1.09
Nonaqueous	13.62	13.15	0.44	~0.3		60.91	-0.35

The data from these calculations show the consistency of this set for percent detected. The most significant observation is the value for the buffering capacity. For HAN analysis, both methods use ~0.5 grams and are similar. But, for the titration of nitric acid, there is a very large difference in buffering capacity which would make the aqueous method more susceptible to error. This larger buffering capacity is the function of the larger sample size, 20 to 30 grams versus the smaller sample for nonaqueous method (0.5 g) .

Based on the optimized nonaqueous titrimetry technique, ionic strengths can be calculated to determine conformity with accepted practices. The molarities and therefore ionic strengths which were used for this study for liquid propellant, LP 1846, containing HAN at 8.94 M and TEAN at 1.34 M are as follows:

$$\begin{aligned}
 &\text{Using 0.6 gram of LP or 0.42 ml at a density of 1.42,} \\
 &\text{HAN} = \text{ml [M/ (1000 ml/L) = 0.00375 moles} \\
 &\text{TEAN} = 0.00056 \text{ moles}
 \end{aligned} \tag{5}$$

$$\begin{aligned}
 &\text{Diluted in 50 ml,} \\
 &\text{HAN} = [ (\text{moles} \times 1000 \text{ ml/L) } / \text{ml}] = 0.074 \text{ M} \\
 &\text{TEAN} = 0.011 \text{ M}
 \end{aligned} \tag{6}$$

$$\begin{aligned}
 &\text{Ionic strength, } \mu, \\
 &\mu = 1/2 \sum \mu z^2 = 0.085
 \end{aligned} \tag{7}$$

In the concentrated systems, the curves are more asymmetrical (figs. 2 through 10) which are a function of the buffering capacity maximum,  $b$ , being close to the titration end-point. In the symmetrical curves from the dilute systems, the maximums are far removed from the end-point.

Another illustration of buffering as expressed by Kolthoff and Sandell states that in any weakly acidic solution,  $HA$ , the equilibrium is determined by the magnitude of the ionization constant:

$$[H^+][A^-] / [HA] = K_a \quad (8)$$

and

$$[H^+] = \{ [HA] / [A^-] \} \times K_a$$

If a mixture of a weak acid and its salt ( $KA$ ) is considered, it is a strong electrolyte. The concentration of  $A^-$  ions sent into solution by the salt is practically equal to the molecular concentration of the salt which represses the dissociation of the acid,  $HA$ , since it furnishes the common ion  $A^-$ . Frequently, this repression is so great that all of the acid present can be considered to be in the undissociated form (ref 8).

During the course of this investigation, it was observed that the break in the titration curve for AN, EAN, and DEAN, which are also weak acids, are not resolved from TEAN. Since these contaminants could have a large impact on the TEAN determination, an alternative method for TEAN is necessary. Due to the observed limitations of titrimetry, ion and supercritical fluid chromatography are being investigated for the analysis of HAN, TEAN, EAN, DEAN, and AN. This does not preclude the use of titrimetry. Potentiometric titrimetry appears to provide accurate HAN analysis and an indication of other contaminants in the titration break for TEAN when a first derivative output is used. The use of titrimetry will be used until such time that either IC or SFC or both are optimized. Both IC and SFC instrument manufacturers \* were consulted on the above-mentioned analyses which has resulted in feasible approaches for further investigation. One of the important considerations which led to SFC was environmental impact and waste disposal.

---

\* J. Krol, Waters Division of Millipore, 1988  
M. Ashraf-Khorassani of Suprex Corporation, 1989

A capability was developed on IC to separate HAN, TEAN, DEAN, EAN, and AN. Current investigation is in progress to determine the precision and accuracy of this technique. Preliminary results indicate that IC is a viable method. An IC-Pak TM column and conductivity detector were used for the separation and detection. A 20 to 100  $\mu$ l size sample of diluted liquid propellant (3/5000 in water) was injected using 1 to 4 mm  $\text{HNO}_3$ /0-15% methanol as the eluant. The separation using this technique resolves all the above-mentioned species (refs 2 and 3). The nitric acid molarity effects the speed of elution of the HAN and AN primarily and this methanol concentration, primarily the TEAN. This combination can be used to vary the elution times of the components of interests for a particular application.

As a result of recent progress in the field of SFC, this chromatographic procedure was also reviewed. Conformity to stricter environmental and waste disposal regulations have required a search for compatible techniques which would permit achievement of analytical goals and compliance with these new directives. This method has been shown to provide the desired results and merits further investigation. The chromatographic column used in this separation for all the species of interest was 10 cm x 1.0 mm x 5  $\mu$ m methyl deltabond with supercritical carbon dioxide plus 0.3 formic acid as eluant. The sample of LP 1846 was extracted on celite with supercritical carbon dioxide plus 0.3 formic acid prior to analysis (ref 10).

## CONCLUSIONS

The final selection of analytical techniques for this program is still under investigation. The uses of titrimetry have been illustrated as well as the advantages of chromatography. The development of the later has far-reaching utility and applicability to many facets of this program. In conclusion, all of these techniques have their advantages and disadvantages and must be weighed for the particular application.

The review of analytical methodologies which are applicable for monitoring the storage of liquid propellants has led to the investigation of many diverse techniques. Any or all of these have applicability for the particular requirements. This investigation has endeavored to provide some incite into their advantages and disadvantages.

Table 1. Determination of nitric acid in spiked propellant samples by nonaqueous titration

<u>Sample</u>	<u>%Acid added</u>	<u>%Acid present</u>	<u>Total % acid</u>	<u>% Acid found</u>
Ethanol	0.03	0.0	0.03	0.02
	0.03		0.03	0.02
	0.05		0.05	0.04
	0.05		0.05	0.04
	0.13		0.13	0.12
	0.13		0.13	0.12
LP 1846 No.1	0.0	0.08	0.08	0.09
	0.0			0.07
	0.04		0.12	0.10
	0.04		0.12	0.11
	0.05		0.13	0.14
	0.06		0.14	0.17
	0.14		0.22	0.21
	0.16		0.24	0.22
	0.33		0.41	0.43
	0.33		0.41	0.41
	0.0		0.11	0.11
	0.0			0.11
LP 1846 No. 2	0.06	0.11	0.17	0.15
	0.16		0.27	0.23
	0.32		0.43	0.41
	0.32		0.43	0.41
LP 1846 No. 3	0.0	0.10	0.10	0.11
				0.09
				0.11
				0.11
	0.03		0.13	0.13
				0.13
				0.13
				0.12
	0.07		0.17	0.17
				0.16
				0.17
				0.16
LP 1846 No. 4	0.03	0.11	0.14	0.00*

Used 2.0 g sample resulting in complete loss in recovery; all other samples are 0.5 g.

Table 2. Determination of nitric acid in spiked propellant samples by aqueous titrations

<u>Sample</u>	<u>Sample wgt,g</u>	<u>% Acid added</u>	<u>%Acid present</u>	<u>% Nitric acid actual total</u>	<u>%Acid found</u>
DI H2O	1.0	0.0	0.00	0.00	0.01
	1.0	0.00	0.00	0.00	0.01
	1.0	0.065	0.065	0.065	0.069
LP 1846 No.1	0.50	0.00	--	--	0.06
	0.47		--	--	0.07
	0.49		--	--	0.06
	0.99		--	--	0.02
	1.03				a,b
	1.97		--	--	0.02
	4.27		--	--	0.04
	4.31		--	--	0.03
	4.29		--	--	0.04
LP 1846 No. 2	0.48	0.03	0.04	0.07	0.05
	0.49	0.03		0.07	0.07
	0.49	0.03		0.07	0.06
	0.50	0.065	0.04	0.10	0.09
	0.50	0.065	0.04	0.10	0.09
	0.49	0.065		0.10	0.11
	0.48	0.17	0.04	0.21	0.22
	0.48	0.29	0.04	0.33	0.31
	0.48	0.28		0.32	0.32
	0.49	0.33		0.37	0.35
LP 1846 No. 3	1.02	0.06	0.04	0.10	a,b
LP 1846 No. 4	4.26	0.004	0.04	0.04	b
	4.26	0.004		0.04	b
	4.23	0.03	0.04	0.07	b
LP 1846 No. 5	28.65 <sup>c</sup>	0.00	0.04	0.04	a,b
	28.53 <sup>c</sup>	0.03		0.07	a,b
	28.61 <sup>c</sup>	0.00		0.04	0.03
	28.53 <sup>c</sup>	0.03		0.07	0.06

<sup>a</sup>Nondiscernible end points using large samples.

<sup>b</sup>Determinations using S- curves.

<sup>c</sup>Total volume is 60 ml; all others are 50 ml.

Table 3. Effect of diluent blank in nitric acid study

<u>Diluent</u>	<u>Titration</u>	<u>Nitric acid</u>	
		<u>Present</u>	<u>Found</u>
Methanol HPLC grade	Nonaqueous	0.00	0.01
Ethanol 1 200 proof	Nonaqueous	0.00	0.01
Ethanol 2 200 proof	Nonaqueous	0.00	0.12
DI water 1	Aqueous	0.00	0.01
DI water 2	Aqueous	0.00	0.03

Table 4. Comparison of aqueous versus nonaqueous titrations for HAN and TEAN

<u>Sample</u>	<u>Aqueous (NaOH)</u>		<u>Nonaqueous (TBAH)</u>	
	<u>% HAN</u>	<u>% TEAN</u>	<u>HAN</u>	<u>% TEAN</u>
LP-2	58.69	21.61	60.42	20.47
	58.67	21.00	60.45	20.42
	58.64	21.48	60.39	20.43
	58.61	21.46		
	58.65	21.76		
	58.65 $\pm$ 0.02	21.46 $\pm$ 0.19	60.42 $\pm$ 0.02	20.44 $\pm$ 0.02
LP-1846-01 (LP 3)	58.76	21.49	60.54	20.90
	58.85	21.51	60.60	20.95
	58.81	21.76	60.64	20.82
	58.71	22.18		
	58.84	22.35		
	58.70	21.75		
	58.78 $\pm$ 0.06	22.01 $\pm$ 0.33	60.59 $\pm$ 0.04	20.89 $\pm$ 0.05

## NOTE:

(1) All values are corrected for free nitric acid and ammonium nitrate.

(2) LP-2 is lot ABY87FS2C013.



Table 5. Comparison of aqueous and non aqueous synthetic  
1846 HACL/TEACL titrations

<u>Titration type</u>	<u>Diluent</u>	<u>Percent acetone</u>	<u>Reaction time, min</u>	<u>Titrant HA · CL</u>	<u>Volume, mls TEA · CL</u>
Nonaqueous	Ethanol	1.0	0.0	14.91	17.28
				<u>14.87</u>	<u>17.23</u>
				<u>+0.02</u>	<u>+0.025</u>
		5.0	0.0	14.95	17.33
				<u>14.93</u>	<u>17.31</u>
				<u>+0.03</u>	<u>+0.01</u>
		10.0	0.0	14.93	17.32
				<u>14.87</u>	<u>17.26</u>
				<u>+0.03</u>	<u>+0.03</u>
		10.0	15.0	<u>14.89</u>	<u>17.26</u>
	Actual concentration		60.8% HAN	19.2%	TEAN
	Experimental		61.4% HAN	21.6%	TEAN
Aqueous Water		1.0	0.0	17.19	20.13
				<u>17.14</u>	<u>20.06</u>
				<u>+0.025</u>	<u>+0.035</u>
		5.0	0.0	17.11	20.03
				<u>17.20</u>	<u>20.10</u>
				<u>+0.045</u>	<u>+0.035</u>
		10.0	0.0	17.14	20.03
				<u>17.15</u>	<u>20.02</u>
				<u>+0.005</u>	<u>+0.005</u>
		10.0	15.0	<u>17.14</u>	<u>20.04</u>
	Actual concentration		60.8% HAN	19.2%	TEAN
	Experimental		58.5% HAN	21.9%	TEAN

NOTE:

High purity hydroxylammonium hydrochloride and triethanolammonium hydrochloride was used to prepare solutions containing the same amount of cations (hydroxylamine or triethanolamine) which would be present in the nitrate solutions of LP 1846.

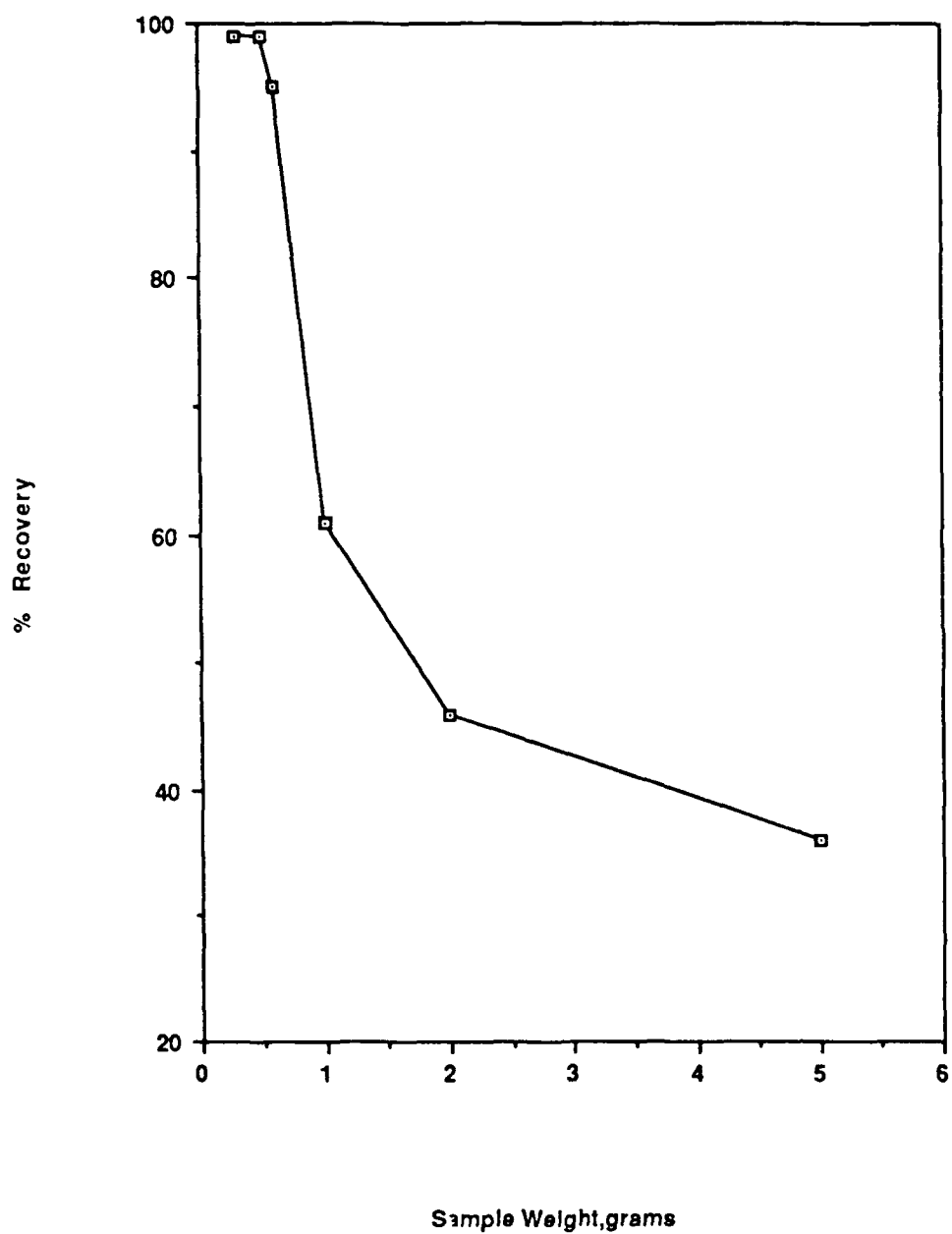


Figure 1. Nitric acid titration, % recovery versus sample weight

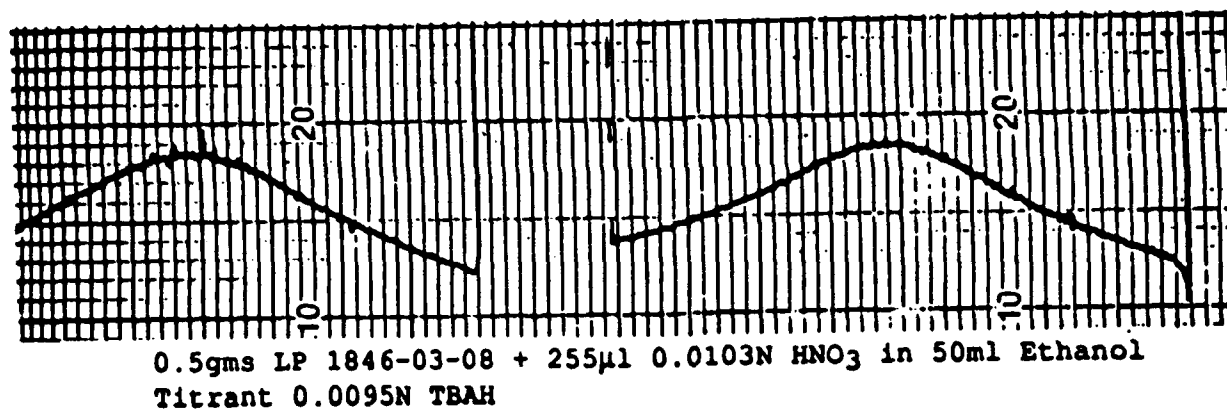


Figure 2. Nonaqueous titration derivative curves with optimized sample size

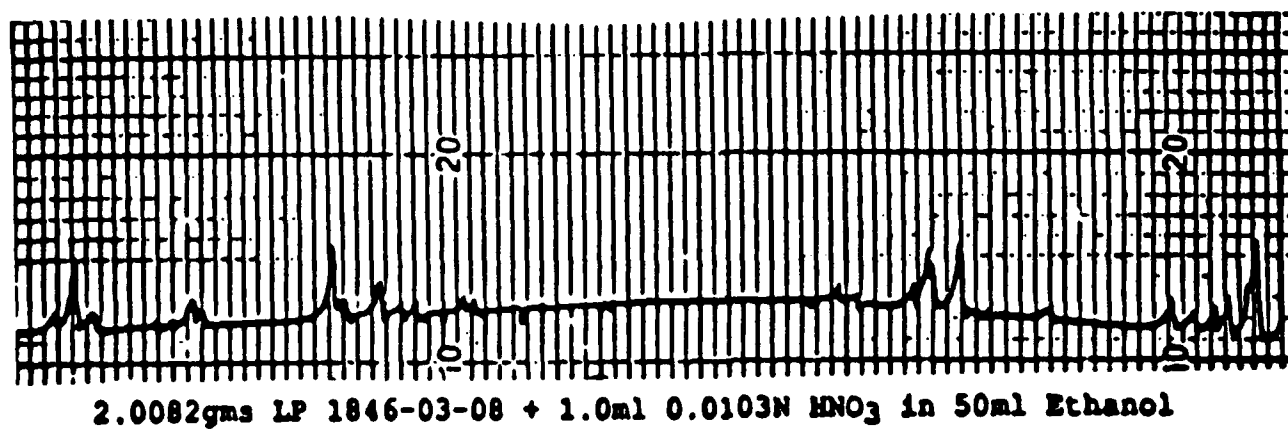


Figure 3. Nonaqueous titration derivative curve with concentrated sample

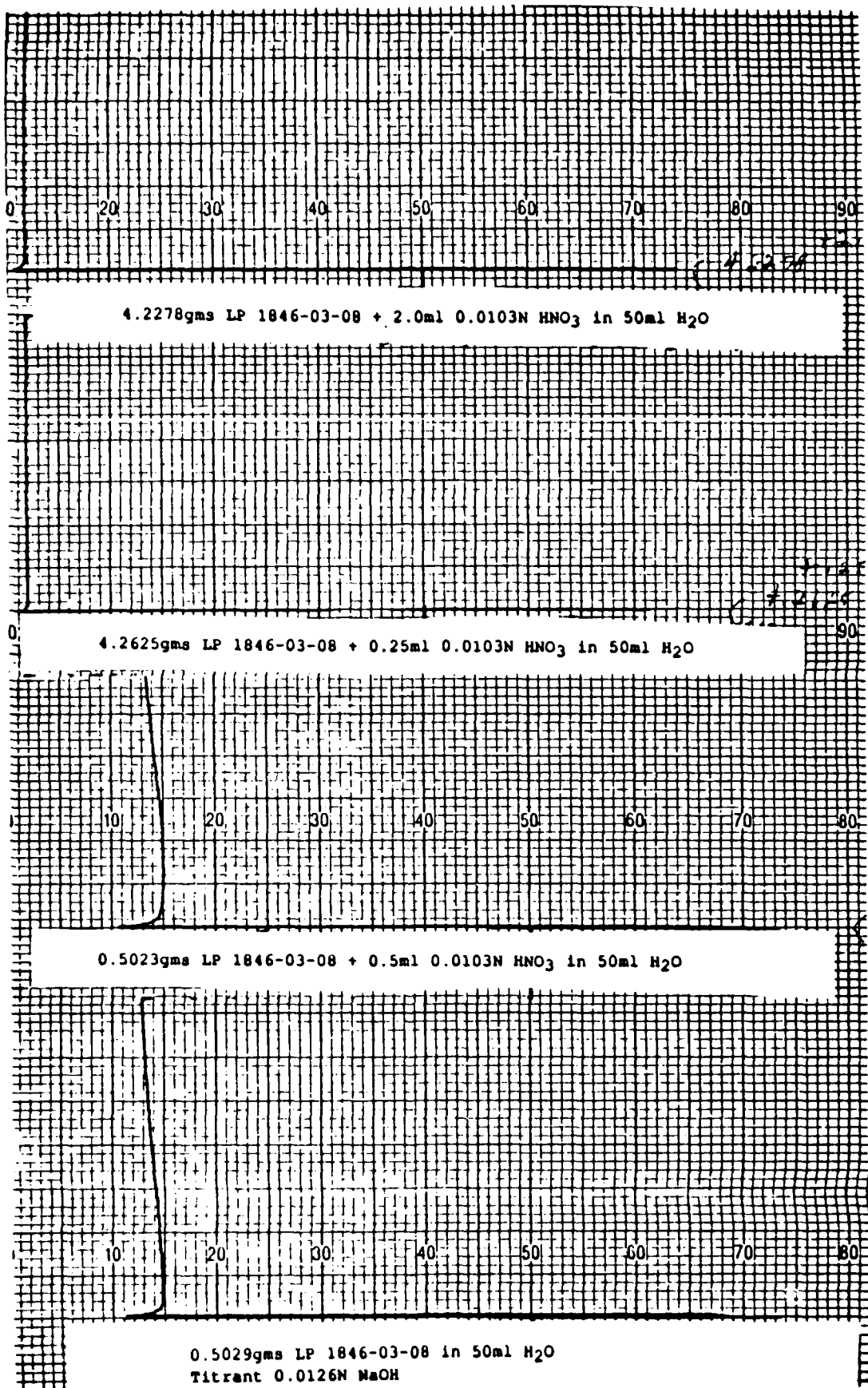


Figure 4. Effect of sample concentration in aqueous titration derivative curves

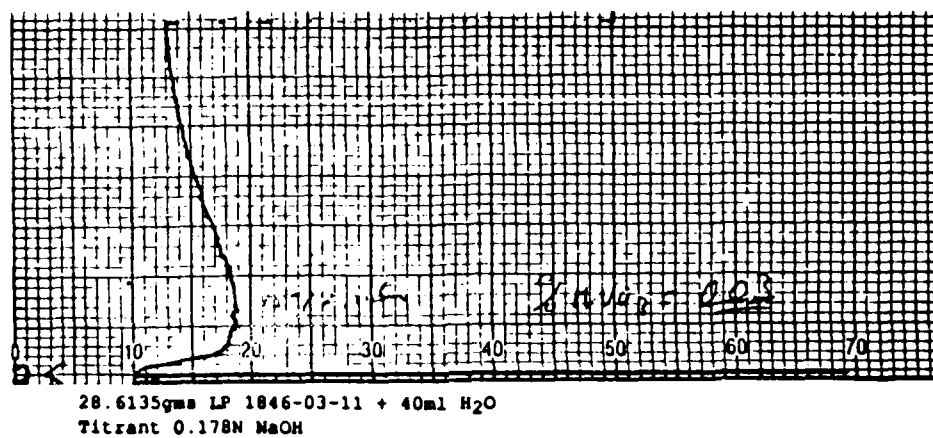
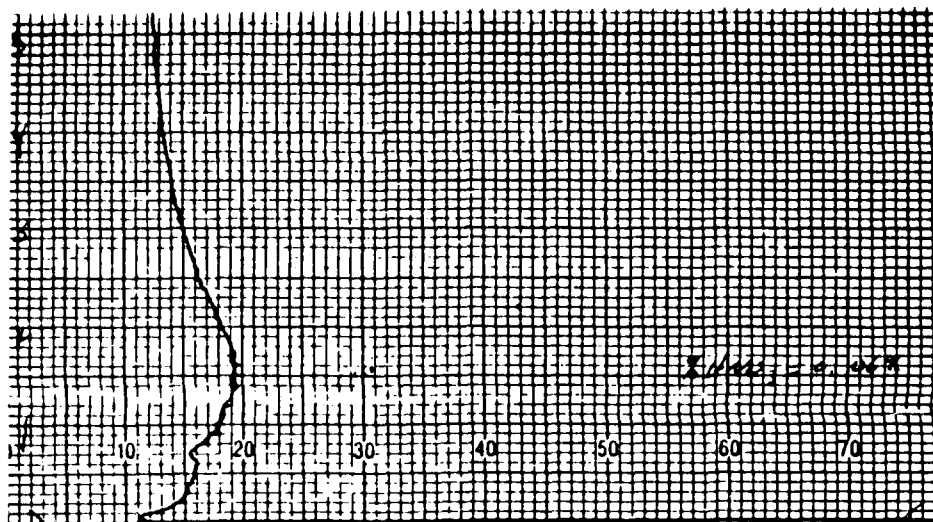


Figure 5. Aqueous titration derivative curves

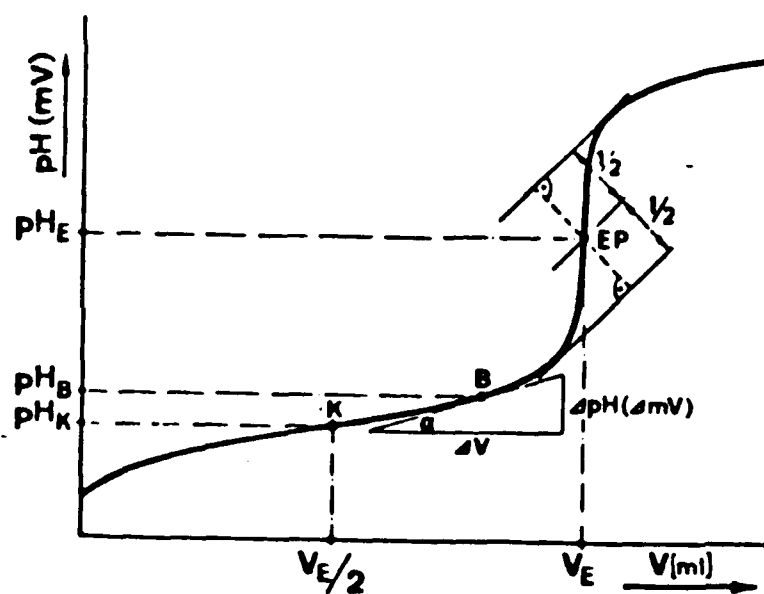


Figure 6. Potentiometric titrations S-curve

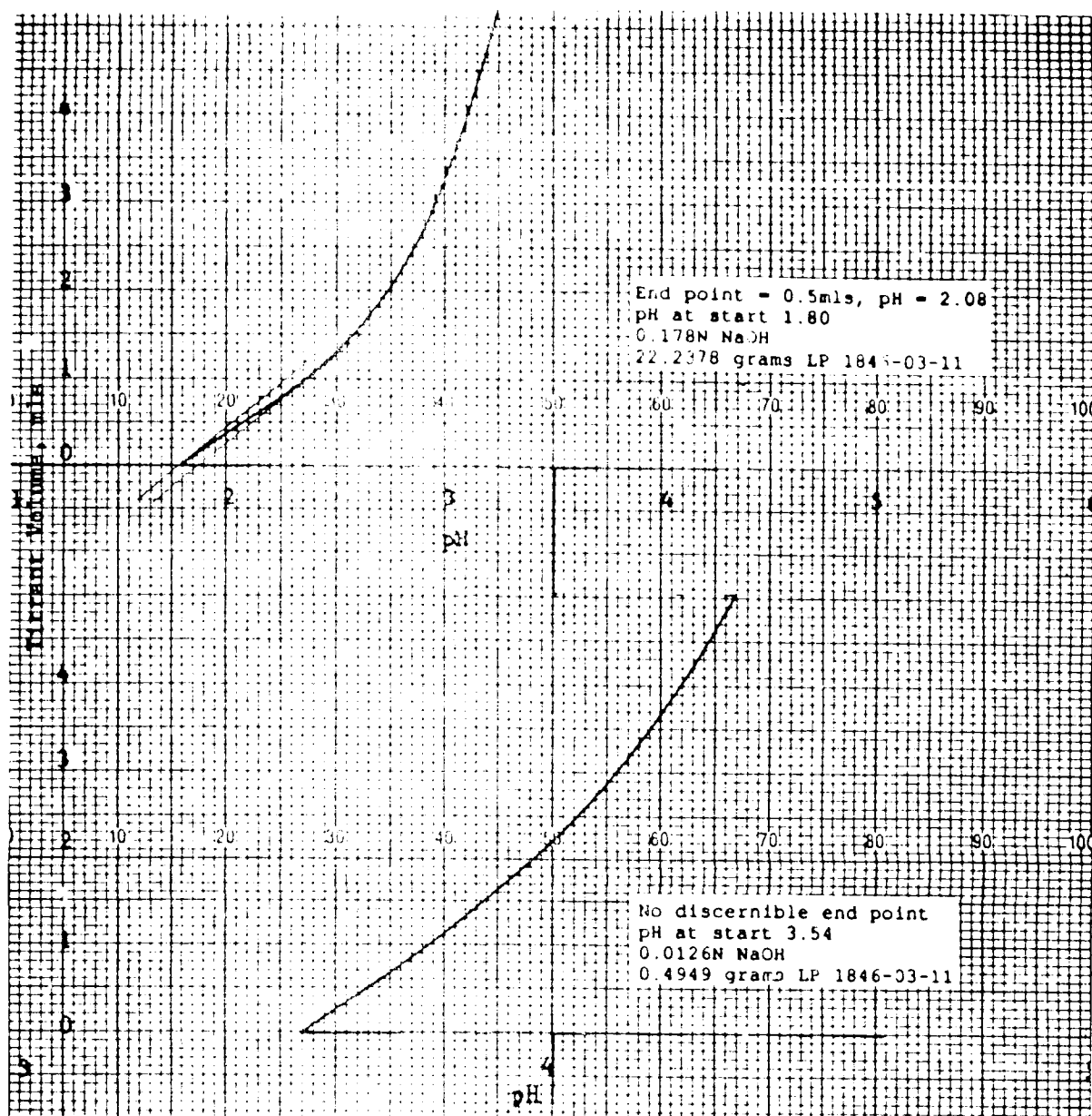


Figure 7. Aqueous titration S-curves for nitric acid

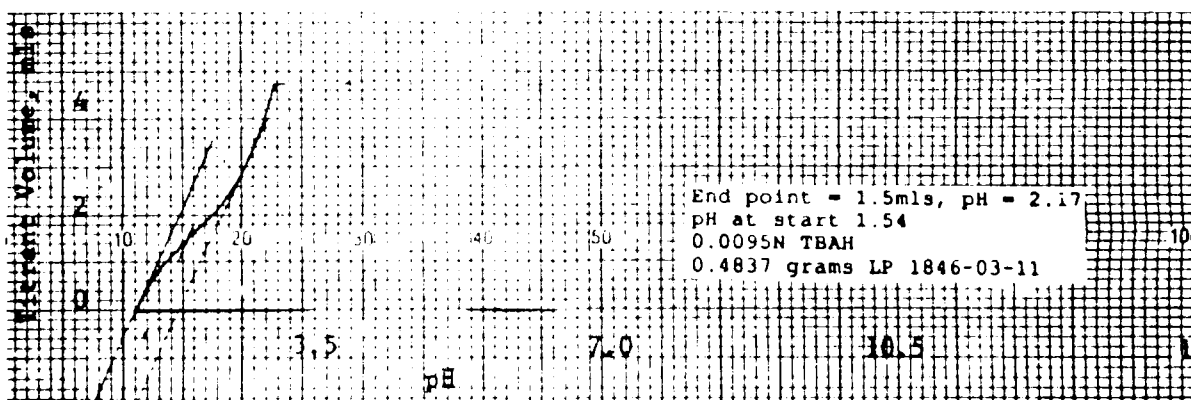


Figure 8. Nonaqueous titration S-curves for nitric acid

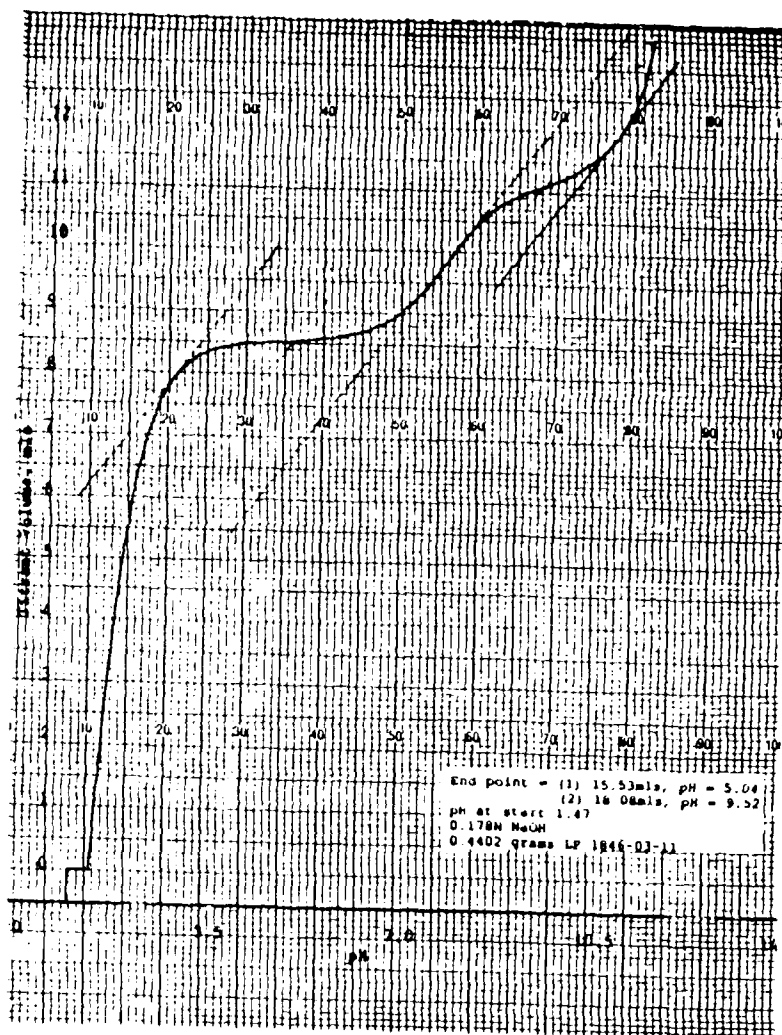


Figure 9. Aqueous titration S-curves for HAN/TEAN

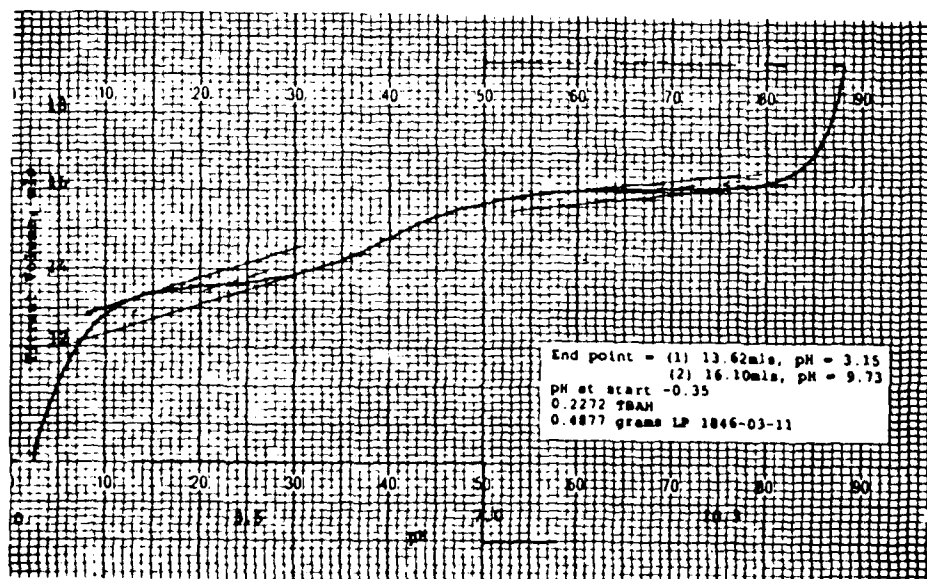


Figure 10. Nonaqueous titration S-curve HAN/TEAN

## REFERENCES

1. Seals, W.O., "Vulnerability Studies for a Liquid Propellant," Technical Report BRL-TR-2646, 1984.
2. Griff, S.P., Seals, W.O., and Doyle, G., "The Effect of the Temperature and Transition Metals on the Stability of LP 1846," 4th International Gun Propellant Symposium, ARDEC, Picatinny Arsenal, NJ., November 1988, Submitted for CPIA Publication.
3. Griff, S.P., Doyle, G., and Seals W.O., "Stability of Liquid Propellants Under Long-Term Storage," JANNAP Propulsion Systems Hazards Subcommittee Meeting, Brooks AFB, San Antonio, Texas, Febuary 1989, Contractor report ARAED-CR-89002, Picatinny Arsenal, NJ, April 1989.
4. Decker, M., Klien, N., and Freedman, E., " Titrimetric Analysis of Liquid Monopropellants, " JANNAP Propellant Characterization Subcommittee Meeting, 1986.
5. Decker, M. and Freedman, E., " Analysis of Liquid Monopropellants by Titration, " JANNAP 23rd Combustion Meeting, 1986.
6. Kolthoff, I.M. and Elving, P.J., "Treatise on Analytical Chemistry," Part I, Volume I, Chapters 10, 11, 12, and 13, Interscience, New york, 1959.
7. Bates, R.G., Determinations of pH, Theory and Practice, Wiley, New York, 1964.
8. Kolthoff, I.M. and Sandell, E.B., Textbook of Quantitative Inorganic Analysis, 3rd Edition, Macmillan, New York, 1952.



## DISTRIBUTION LIST

Administrator  
Defense Technical Info Center  
ATTN: Accessions Division (12)  
Cameron Station  
Alexandria, VA 22304-6145

HQDA (SARD-TR)  
Washington, D.C. 20310-0001

Commander  
U.S. Army Materiel Command  
ATTN: AMCDRA-ST  
5001 Eisenhower Avenue  
Alexandria, VA 22333-0001

Commander  
U.S. Army Laboratory Command  
ATTN: AMSLC-DL  
Adelphi, MD 20783-1145

Director  
Benet Weapons Laboratory  
Armament RD&E Center  
U.S. Army AMCCOM  
ATTN: SMCAR-CCB-TL  
Watervliet, NY 12189-4050

Commander  
U.S. Army Armament, Munitions and Chemical Command  
ATTN: SMCAR-ESP-L  
Rock Island, IL 61299-5000

Commander  
U.S. Army Aviation Systems Command  
ATTN: AMSAV-DACL  
4300 Goodfellow Blvd.  
St. Louis, MO 63120-1798

Director  
U.S. Army Aviation Research and Technology Activity  
Ames Research Center  
Moffett Field, CA 94035-1099

Commander  
U.S. Army Missile Command  
ATTN: AMSMI-RD-CS-R (DOC)  
Redstone Arsenal, AL 35898-5010

Commander  
U.S. Army Tank-Automotive Command  
ATTN: AMSTA-TSL (Technical Library)  
Warren, MI 48397-5000

Director  
U.S. Army TRADOC Analysis Command  
ATTN: ATAA-SL  
White Sands Missile Range, NM 88002-5502

Commandant  
U.S. Army Infantry School  
ATTN: ATSH-CD-CSO-OR  
Fort Benning, GA 31905-5660

Air Force Armament Laboratory  
ATTN: AFATL/DLODL  
Eglin AFB, FL 32542-5000

Director  
USAMSAA  
ATTN: AMXSY-D  
AMXSY-MP, H. Cohen  
Aberdeen Proving Ground, MD 21005-5066

Commander  
USATECOM  
ATTN: AMSTE-TD  
Aberdeen Proving Ground, MD 21005-5066

Commander  
CRDEC  
AMCCOM  
ATTN: SMCCR-RSP-A  
Aberdeen Proving Ground, MD 21010-5423

Commander  
CRDEC  
AMCCOM  
ATTN: SMCCR-MSI  
SMCCR-MU  
Aberdeen Proving Ground, MD 21010-5423

Director  
VLAMO  
ATTN: AMSLC-VL-D  
Aberdeen Proving Ground, MD 21005-5066

Director  
Defense Advanced Research Projects Agency  
ATTN: J. Lupo  
J. Richardson  
1400 Wilson Boulevard  
Arlington, VA 22209

HQDA  
ATTN: SARD-ZT, G. Singley  
SARD-TT, I. Szkrybalo  
SARD-TC, C. Church  
D. Zimmerman  
Washington, D.C. 20310

HQ  
U.S. Army Materiel Command  
ATTN: AMCICP-AD, B. Dunetz  
5001 Eisenhower Avenue  
Alexandria, VA 22333-0001

Commander  
Materials Technology Laboratory  
U.S. Army Laboratory Command  
ATTN: SLCMT-MCM-SB, M. Levy  
Watertown, MA 02172-0001

Commander  
USACECOM R&D Technical Library  
ATTN: ASQNC-ELC-I-T, Myer Center  
Ft. Monmouth, NJ 07703-5000

Commander

Armament RD&E Center

U.S. Army AMCCOM

ATTN: SMCAR-IMI-I (5)

SMCAR-AEE-B, D. Downs

SMCAR-AEE-BR, B. Brodman

W. Seals

A. Beardell

SMCAR-AEE-W, N. Slagg

SMCAR-AEE, A. Bracuti

J. Lannon

M. Gupta

J. Salo

D. Chieu

SMCAR-FSS-D, L. Frauen

SMCAR-FSA-S, H. Liberman

SMCAR-FSS-DA, J. Feneck

R. Kopmann

J. Irizarry

SMCAR-CCS-C, T. Hung

SMCAR-SFS

Picatinny Arsenal, NJ 07806-5000

Director

Benet Weapons Laboratory

Armament RD&E Center

U.S. Army AMCCOM

ATTN: SMCAR-CCB, L. Johnson

SMCAR-CCB-S, F. Heiser

SMCAR-CCB-DS, E. Conroy

A. Graham

Watervliet, NY 12189-4050

Commander

U.S. Army Harry Diamond Laboratories

ATTN: SLCHD-TA-L

2800 Powder Mill Road

Adelphi, MD 20783-1145

Commander

U.S. Army Belvoir Research and Development Center

ATTN: STRBE-WC

Technical Library (Vault) B-315

Fort Belvoir, VA 22060-5606

Commander  
U.S. Army Research Office  
ATTN: Technical Library  
P.O. Box 12211  
Research Triangle Park, NC 27709-2211

Commandant  
U.S. Army Field Artillery School  
ATTN: ATSF-CMW  
ATSF-TSM-CN, J. Spicer  
Fort Sill, OK 73503

Commandant  
U.S. Army Armor Center  
ATTN: ATSB-CD-MLD  
Fort Knox, KY 40121

Commander  
Naval Surface Weapons Center  
ATTN: D.A. Wilson (Code G31)  
Dahlgren, VA 22448-5000

Commander  
Naval Surface Weapons Center  
ATTN: J. East (Code G33)  
Dahlgren, VA 22448-5000

Commander  
U.S. Naval Surface Weapons Center  
ATTN: O. Dengrel  
K. Thorsted  
Silver Spring, MD 20902-5000  
Commander (Code 3247)  
Naval Weapons Center  
Gun Systems Branch  
China Lake, CA 93555-6000

Superintendent  
Naval Postgraduate School  
Department of Mechanical Engineering  
ATTN: Code 1424, Library  
Monterey, CA 93943

AFOSR/NA (L. Caveny)  
Building 410  
Bolling AFB  
Washington, DC 20332

Commandant  
USAFAS  
ATTN: ATSF-TSM-CN  
Fort Sill, OK 73503-5600

Director  
Jet Propulsion Laboratory  
ATTN: Technical Library  
4800 Oak Grove Drive  
Pasadena, CA 91109

Director  
National Aeronautics and Space Administration  
ATTN: MS-603, Technical Library  
MS-86, Dr. Povinelli  
21000 Brookpark Road  
Lewis Research Center  
Cleveland, OH 44135

Director  
National Aeronautics and Space Administration  
Manned Spacecraft Center  
Houston, TX 77058

Central Intelligence Agency  
Office of Central Reference  
Dissemination Branch  
Room GE-47 HQS  
Washington, DC 20502

Central Intelligence Agency  
ATTN: Joseph E. Backofen  
HQ Room 5F22  
Washington, DC 20502

Calspan Corporation  
ATTN: Technical Library  
P.O. Box 400  
Buffalo, NY 14225

General Electric Ordnance System Division

ATTN: J. Mandzy, OP43-220

R.E. Mayer

H. West

W. Pasko

R. Pate

I. Magoon

J. Scudiere

Minh Luu

100 Plastics Avenue

Pittsfield, MA 01201-3698

General Electric Company

Armament Systems Department

ATTN: D. Maher

Burlington, VT 05401

Honeywell, Inc.

ATTN: R.E. Tompkins

MN38-330

10400 Yellow Circle Drive

Minnetonka, MN 55343

IITRI

ATTN: Library

10 West 35th Street

Chicago, IL 60616

Olin Chemicals Research

ATTN: David Gavin

P.O. Box 586

Cheshire, CT 06410-0586

Paul Gough Associates

ATTN: Paul Gough

1048 South Street

Portsmouth, NH 03801-5423

Safety Consulting Engineer

ATTN: Mr. C. James Dahn

5240 Pearl Street

Rosemont, IL 60018

Sandia National Laboratories  
ATTN: R. R. Rychnovsky, Div 8152  
S. Griffiths, Div 8244  
R. Carling, Div 8152  
P.O. Box 969  
Livermore, CA 94551-0969

Science Applications, Inc.  
ATTN: R. Edelman  
23146 Cumorah Crest  
Woodland Hills, CA 91364

Science Applications International Corporation  
ATTN: Dr. F.T. Phillips  
Dr. Fred Su  
10210 Campus Point Drive  
San Diego, CA 92121

Science Applications International Corporation  
ATTN: Norman Banks  
4900 Waters Edge Drive  
Suite 255  
Raleigh, NC 27606

Sundstrand Aviation Operations  
ATTN: Mr. Owen Briles  
P.O. Box 7202  
Rockford, IL 61125

Veritay Technology, Inc.  
ATTN: E.B. Fisher  
4845 Millersport Highway  
P.O. Box 305  
East Amherst, NY 14051-0305

Director  
The Johns Hopkins University  
Applied Physics Laboratory  
Johns Hopkins Road  
Laurel, MD 20707



Director  
CPIA  
The Johns Hopkins University  
ATTN: T. Christian  
Technical Library  
Johns Hopkins Road  
Laurel, MD 20707

University of Illinois at Chicago  
ATTN: Professor Sohail Murad  
Department of Chemical Engineering  
Box 4348  
Chicago, IL 60680

University of Maryland at College Park  
ATTN: Professor Franz Kasler  
Department of Chemistry  
College Park, MD 20742

University of Missouri at Columbia  
ATTN: Professor R. Thompson  
Department of Chemistry  
Columbia, MO 65211

University of Michigan  
ATTN: Professor Gerard M. Faeth  
Department of Aerospace Engineering  
Ann Arbor, MI 48109-3796

Director  
Ballistic Research Laboratory  
ATTN: AMXBR-OD-ST  
Aberdeen Proving Ground, MD 21005-5423

Commander  
U.S. Army Armament, Munitions and Chemical Command  
ATTN: AMSMC-GCL (D)  
Picatinny Arsenal, NJ 07806-5000

University of Missouri at Columbia  
ATTN: Professor F.K. Ross Research Reactor  
Columbia, MO 65211

University of Missouri at Kansas City  
Department of Physics  
ATTN: Professor R.D. Murphy  
1110 East 48th Street  
Kansas City, MO 64110-2499

Pennsylvania State University  
Department of Mechanical Engineering  
ATTN: Professor K. Kuo  
University Park, PA 16802

Princeton Combustion Research Laboratories, Inc.  
ATTN: N.A. Messina  
M. Summerfield  
4275 U.S. Highway One North  
Monmouth Junction, NJ 08852

University of Arkansas  
Department of Chemical Engineering  
ATTN J. Havens  
227 Engineering Building  
Fayetteville, AR 72701

University of Delaware  
Department of Chemistry  
ATTN: MR. James Cronin  
Professor Thomas Brill  
MR. Peter Spohn  
Newark, DE 19711

University of Texas at Austin  
Bureau of Engineering Research  
ATTN: BRC EME133, Room 1.100  
H. Fair  
10100 Burnet Road  
Austin, TX 78758

Olin Corporation  
ATTN: Victor A. Corso  
DR. Ronald L Dotson  
24 Science Park  
New Haven, CT 06511



DEPARTMENT OF THE ARMY  
U.S. ARMY ARMAMENT RESEARCH, DEVELOPMENT AND ENGINEERING CENTER  
PICATINNY ARSENAL, NEW JERSEY 07806-5000

# ERRATA

REPLY TO  
ATTENTION OF

SMCAR-IMI-I

11 September 1990

SUBJECT: Report Corrections

Defense Technical Information Center  
Selection Section  
DTIC-FDA-C  
Cameron Station, VA 22304-6145

Corrections should be made to the following reports:

1. ADA 226186 - author's name should be S. P. Griff not J. P. Griff
2. ADB 147023 - report number should be ARFSD-SP-90002 not  
ARFSD-SP-90004

*Florence Winkelman*

FLORENCE WINKELMAN  
Editorial Services  
STINFO, IMD

AD-A226186